


# *Escali*<sup>®</sup>

**SMARTCONNECT™**

**BODY SCALE**

MODEL # SC200BS

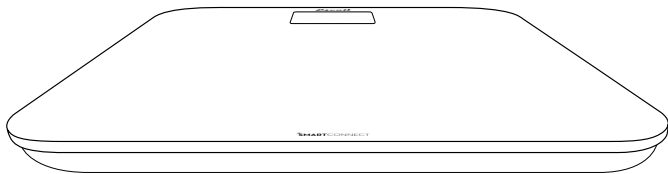

## **OWNER'S MANUAL**

Thank you for your purchase! Please read the owner's manual carefully to ensure proper function of your scale.

**WWW.ESCALI.COM**

# Quick Start Guide

## INITIAL SETUP

1. Insert batteries into the battery compartment located beneath the scale.
2. Download the free Escali SmartConnect app from the Apple App Store. Direct link: [www.escali.com/app](http://www.escali.com/app)

## PAIRING WITH THE ESCALI SMARTCONNECT APP

1. Turn on bluetooth for your iPhone or iPad.
2. Install and open the free Escali SmartConnect app.
3. Create a profile or sign in to your profile.
4. Navigate to the HEALTH screen & select 'WEIGH IN'.
5. Press the button located underneath the body scale to pair devices. Two brackets facing each other “[ ]” on the body scale's display will confirm successful pairing.

### Hardware requirements:

iPhone 4s (or newer), iPad 3rd generation (or newer), iPad Mini (or newer), iPod Touch 5th generation (or newer)

### Software requirements:

iOS 6.0 (or newer)

**Attention!** You must pair the Escali SmartConnect Body Scale with the Escali SmartConnect app once before weight measurements will save in your profile.

## PRODUCT SUPPORT

24/7 Online community support: <https://escali.zendesk.com>

Email support: [feedback@escali.com](mailto:feedback@escali.com)

Phone support: 1-800-467-6408

# Product Guide

## **WEIGHT MEASUREMENTS**

1. Step onto the center of the scale, it will turn on automatically.
2. Step off after weight is displayed. The scale will turn itself off.
3. Press the button underneath the scale to change units.

## **BLUETOOTH® SMART**

Bluetooth® SMART wireless technology allows for the Escali SmartConnect Body Scale to automatically send weight measurements to your iPhone or iPad for enhanced functionality. Download the Escali SmartConnect app from the Apple App Store.

## **INSTANT-ON TECHNOLOGY**

The scale will turn on automatically when weight is applied to the platform.

## **DISPLAY HOLD**

Your weight will stay on display for 15 seconds after the measurement is taken, even if you step off the scale.

## **AUTOMATIC SHUT-OFF**

The scale will turn off after 15 seconds of inactivity to save battery.

## **LARGE BACKLIGHT DISPLAY**

- The backlight will turn on when adding items onto the scale.
- The backlight turns off automatically after 18 seconds of inactivity.

## **PROPER USE & CARE**

- Use and stow the scale at room temperature.
- Use on a flat, level and stable surface.
- The scale is an instrument of precision. Do not submit the scale to shaking, excessive vibration or other rough treatment.
- The display may be affected by electromagnetic disturbances

such as radios, microwaves or cell phones. If such disturbances occur, remove the source of the disturbance & restart the scale.

- Disassembling or tampering with the scale will void the warranty.
- Clean the scale with dish soap and a damp cloth.
- Do not submerge in liquid. The scale is not waterproof.

## WARNING MESSAGES

---- Overload warning.

 Low battery warning.

## PRODUCT SPECIFICATIONS

Model: SC200BS

Capacity: 440 lb. (200 kg.)

Graduation: 0.2 lb. (0.1 kg.)

Dimensions: 12-5/8 inches x 12-5/8 inches x 1 inch

Power Source: 4 AAA batteries (included)

## LIFETIME LIMITED WARRANTY

This scale is guaranteed to be free of defects in materials and workmanship. For more details visit: [www.escali.com/warranties](http://www.escali.com/warranties)

## ESCALI CORPORATION

3203 Corporate Center Drive, Suite 150

Burnsville, MN 55306

**Escali's liability is limited to two (2) times the cost of the product.**

© Escali Corp., Minneapolis, MN, USA. All rights reserved. Escali, the Escali logo and other Escali marks are owned by Escali Corp. and may be registered. The Bluetooth® SMART word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. iPhone, iPad and iPad Mini are trademarks of Apple Inc., registered in the U.S. and other countries.