


LINGLEE

HAIR REMOVAL MACHINE

Model: TFDA28

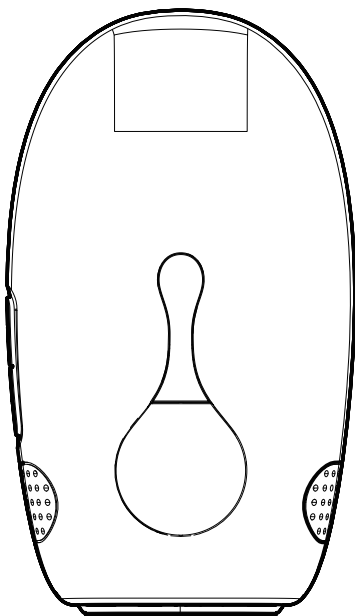

USER MANUAL

Please read the instructions before use and keep for future reference.

Contents

Introduction.....	01
Important Safety Instructions.....	01
Product Description.....	05
Specifications.....	05
Before Use.....	06
Steps for Hair Removal.....	08
Steps for Ice Compress Mode.....	09
Frequency for Reference.....	10
Cleaning and Maintenance.....	11
Trouble Shooting.....	11
Warranty.....	14

Introduction

Hello! Thank you for purchasing a Linglee product.

In order to best use this product and to prevent personal injury, please be sure to read the safety instructions in this manual before use.

Important Safety Instructions

To prevent the risk of burns, fire, electric shock or injury:


- ⊗ Do not use this appliance with wire tangled. Do not twist the power cord during use or storage.
- ⊗ Do not wrap the power cord around the body of the device. Do not apply undue force to the power adapter, socket, or power cord.
- ⊗ If the product is damaged, do not touch any internal parts.
- ⊗ Do not use this product around the eyes.
- ⊗ We recommend wearing protective glasses for every use. Although the light will not directly affect your eyes, the high-frequency flashing light may affect your vision.
- ⊗ When you need back hair removal or ice compress, you may need an assistant or partner. It may be risky if you use it alone.
- ⊗ If your hair is longer than 1mm, please shave before use.
NOTE: Excessively long hair absorbing heat on the surface of the skin will cause skin redness and swelling, and scorched hair will produce peculiar smells, which will affect your experience of using the product.
- ⊗ Prevent children under 15 years of age from using this product.
- ⊗ If the power cord is damaged please contact consumer support.
- ⊗ **The following people should not use the product.**
 - If you are pregnant, menstruating or lactating.
 - If you have body surface scars or diabetes.
 - If you have sensitive skin or light-sensitive skin.
 - If you have a severe heart disease or hypertension.

- If you have malignant tumor or folliculitis.
 - If you have infectious skin disease or burned skin.
 - ⊗ **Do not use it including but not limited to the following conditions.**
 - If you have abnormal collagen metabolism, including a history of scars or poor wound healing.
 - If you have vascular diseases, such as varicose veins at the irradiation site.
 - If the skin is sensitive to light and is prone to rash or allergies.
 - If the irradiated site is infected, has rash, skin lesions, hematoma, etc.
 - If you have suppressed immune response diseases (including HIV / AIDS infection etc.)
 - If you have hemorrhagic disease.
 - If you have used chemical methods for hair removal within the past 30 days.
 - If you have used skin whitening products or vitamins within the past 30 days.
 - Please be cautious when using other cosmetics with Vitamin A
 - Cannot be used with essential oil products. Please make sure that the skin on the exposed area is clean and dry with no creams, lubricants, or other cosmetics are used.
- NOTE: If the surface of the skin is not clean and not dry it may affect the use experience of the product, and even damage the product.**
- If you are using drugs such as aspirin.
 - If you have recently used light-sensitive drugs (such as tetracycline, sulfa, etc.) do not use this product for at least 30 days after last use.
 - If you have sunbur, please wait for the skin to repair before using it. Do not use it on sensitive skin.
 - If you have any type of wound that has not healed do not use this device.
 - If you take painkillers it will reduce the skin's sensitivity to heat, and it is recommended to use it carefully before the effect of the medicine is relieved.


- If your skin or hair color is not suitable for this device according to the chart below (P06), do not use this product. (If necessary, please consult your doctor before use)
- ⊗ **Do not use it including but not limited to the following areas.**
 - Inside the oral cavity, around genitals, mucous membranes, etc.
 - Brows and hairline
 - Sites where plastic surgery / surgery was performed
 - Ear
 - Jaw knot
 - Nipple and areola
 - Eyeballs, eyelids and areas around the eyes
 - Inside the belly button
 - Head
 - Large veins
 - Acne area
 - Sites with old tattoos or body painting
 - Spots and moles

To reduce the risk of burns, damage or injury:

- ❗ Do not use antiperspirant after irradiating the underarms.
- ❗ Please confirm and clean the dirt of the light outlet port regularly.
- ❗ Do not put the light outlet port close to items that are easy to demagnetize.
- ❗ Do not use when the vent is blocked.
- ❗ Do not irradiate clothes, hair and other objects.
- ❗ Do not use on pets.
- ❗ Do not use this product around the eyes to avoid eye damage.
- ❗ Do not use it near children. The product will emit strong light, which can easily hurt your eyes.
- ❗ After use, please wipe clean the hair and skin care products that attached to the product.
- ❗ Please set the irradiation intensity according to your skin condition. Start using it from level 1 for the first time.

- ❗ After use, if the skin is slightly hot and red, this is normal. The ice-cool mode can be turned on to relieve skin discomfort.
- ❗ When the back hair removal operation is required, it is recommended to have a partner to help.
- ❗ The product cannot be left unattended after it is turned on. Please turn it off after use.
- ❗ Before use, clean the skin and make sure it is completely dry and free of oily stains.
- ❗ If the using environment temperature differs more than 40°F (from high temperature to low temperature or in reverse), please turn off the device and put it aside for 2 hours before using it.
- ❗ Do not use near flammable and explosive materials (including sprayers, gas stations, etc.)
- ❗ Do not touch the light outlet port after use because the light outlet port generates high temperature.
- ❗ Within 24 hours after irradiation, please avoid direct ultraviolet rays on the irradiated part.
- ❗ If the skin of the irradiated area is not dry or clean it may cause the product to not flash or flash unexpectedly.
- ❗ Do not place the product in the sun or expose it to ultraviolet light for a long time.
- ❗ If the built-in fan does not work (no humming sound is heard) after turning on the power, do not use the product.
- ❗ Do not use it in bathing or showering. 
- ❗ Do not use it near a bathtub, shower tray, washbasin or other utensils containing water. When using this product, please keep the product and the power adapter dry.
- ❗ Warning : Class A equipment is intended for use in an industrial environment. In the documentation for the user, a statement shall be included drawing attention to the fact that there may be potential difficulties in ensuring electromagnetic compatibility in other environments, due to conducted as well as radiated disturbances.

Product Description

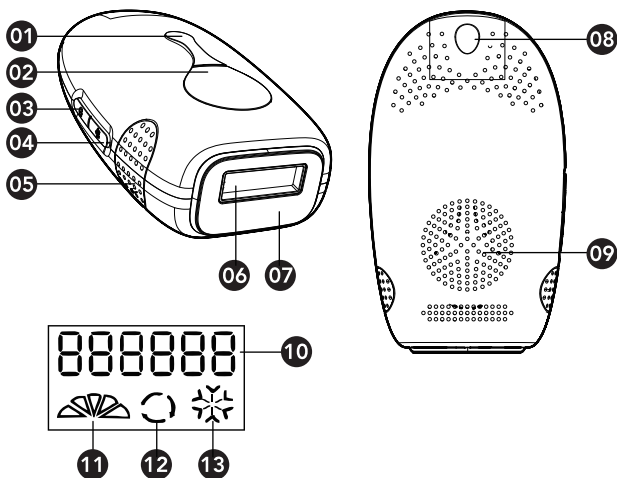

- 01** Manual flash button
- 02** Display screen
- 03** Power/level button
- 04** Ice-cool/automatic button
- 05** Air vent
- 06** Light outlet port
- 07** Ice compress panel
- 08** Charge socket
- 09** Air vent
- 10** Remaining pulse
- 11** Light level indicator
- 12** Automatic flash mode indicator
- 13** Ice-cool mode indicator

Note: Do not block the air vent.

Specifications

Product Name: Hair Removal Machine

Product Mode: TFDA28

Input Power: 12V 3A

Power: 36W

Specifications

Weight: 8.3 oz

Wavelength: 530nm (average)

Operating environment temperature: 41°F-104°F

Recommended environment temperature: 64°F-77°F

Before Use

- Please thoroughly clean the skin and ensure that the skin is hair-free and completely dry and free of any substances.
- Confirm that the light outlet port is free of dirt, hair, etc., and there is no damage or abnormality. If there is dirt, please clean it up to avoid product malfunction, deterioration and skin problems.

Light Energy Levels

- This product is specially developed for home hair removal.
- The product provides 5 levels of light energy with different intensities. You can adjust the level to suit your skin and body hair.
- Please refer to the chart below to determine which level is suitable for your skin and hair color.

Skin Hair	White	Beige	Light brown	Mid brown	Dark brown	Brownish black and darker
White/gray	x	x	x	x	x	x
Red	x	x	x	x	x	x
Light blond	x	x	x	x	x	x
Dark blond / light brown	Level 4/5	Level 4/5	Level 4/5	Level 1/2/3	x	x
Brown	Level 4/5	Level 4/5	Level 4/5	Level 1/2/3	x	x
Dark brown	Level 4/5	Level 4/5	Level 4/5	Level 1/2/3	x	x
Black	Level 4/5	Level 4/5	Level 4/5	Level 1/2/3	x	x

NOTE: x means the product is not suitable nor safe for your skin and / or hair color.

Tips:


- If you feel any pain or discomfort lower the level.
- Compared to white skin, darker skin are more likely to develop skin allergies because the skin absorbs more light. Therefore, people with darker skin are recommended to use a lower level.
- The light on the skin may feel warm after use, but it should never be painful.
- Due to the existence of various reasons, your skin may react differently due to different flash pulses / situations.

Skin Sensitivity Test


If you are using this product for the first time or have recently tanned, you should conduct a skin test on each part to be treated. The skin test is necessary to check whether the skin responds to radiation and to determine the correct light intensity level for each treatment area.




1. Select a part near the area to be treated.


Note: Do not try to use this product on difficult or sensitive parts (foot bare or bone parts).


2. Press and hold the power button  for 2 seconds to turn on, make sure to select level 1.
3. Place the device vertically on the skin so that the light outlet port is in contact with the skin.
4. Press the flash button to release a flash.
5. Slide the device against the skin to the next area to be treated.
6. Increase the setting level by level, release a flash, and then slide the device to the next part. Repeat this operation for all levels within the recommended light intensity range for your skin type.
7. After the skin test, check the skin condition after 2 hours. If there is no abnormality (redness, blisters) in the irradiated area, you can use this product with confidence.


Steps for Hair Removal

1. Before use the product, use a razor to shave the hair in the area that needs to be treated (wet the hair and use the wet shaving method to avoid irritation).
2. After inserting the adapter and connecting the power supply, the display flashes three times and displays "OFF".
3. Find the recommended level in the Light Energy Level chart. If your skin / color hair gets x in the chart, please don't use this device.
4. Wear protective glasses (or sunglasses) during treatment to avoid eye injury.
5. Press and hold the power button  for 2 seconds to turn on the machine.

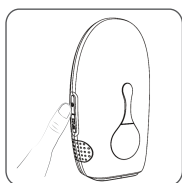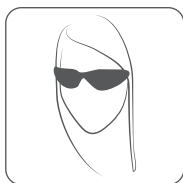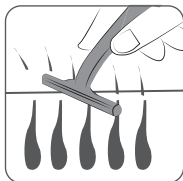

6. At this time, adjust the level according to your skin color, and short press the level button  to adjust the appropriate level. The level you choose will be indicated on the display. It is recommended to start with level 1 for the first use. The higher the level, the better the hair removal effect will be without damaging the skin.
7. Choose to use single flash mode or automatic flash mode.
 - In Single flash mode, after choosing level, press the flash button once and the light will flash once. Place the product vertically on the skin that needs to be treated. If the light outlet port fits the skin completely, the display light will flash. The display light will not flash when it doesn't completely fit.

- In Automatic flash mode, after choosing level, long press the automatic button to start continuous flashes. The cycle icon  will show up on the display. At this time, place the product vertically on the skin that needs to be treated. If the light outlet port fits the skin completely, the display light will flash and stay on.

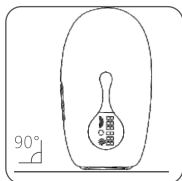**NOTE:**

- There is no ice-cool mode by default (no snow icon on the display). Short press the ice-cool button to start/stop the ice-cool mode. The cooling effect can reduce the burning sensation after use.
- Long press the automatic button to switch between the two modes.
- In Automatic flash mode, the machine will not emit light if it is not in contact with the skin.
- Do not keep flashing on the same part of the skin. Flash once and change to another position to continue removing hair.

Steps for Ice Compress Mode

1. Hold and press the power button for 2 seconds to turn on the machine.
2. Press the ice-cool button to start ice compress mode.
3. Fit the ice compress panel to the part of the hair removal site and move it back and forth to relieve the discomfort of the skin.

NOTE: There's no need to press flash button in ice compress mode.

4. After finishing the treatment you can hold and press the power button for 2 seconds to power off the product or continue to operate in other modes.

NOTE:

- After powering off, due to the high internal temperature, the cooling system will continue to work for a period of time. This is normal. Please wait for it to stop working before unplugging the power supply.

- It is recommended to avoid contact with water within 6 hours after the hair removal operation. After 6 hours, moisturize the treated area.
- Please do not use perfume, facial mask, or cosmetics containing scrub beads after hair removal.
- Within 24 hours after the hair removal operation please avoid UV exposure, and if necessary, please apply sunscreen products. If you feel unwell, please consult a doctor.
- Sunbathing is likely to cause skin darkening or sunburn. Irradiation at this time will cause redness, swelling, inflammation, etc. Therefore, during the use of this product, do not expose your skin to the sun. Wait at least 1 week after your sunburn has recovered before using this device.
- It is strictly forbidden to expose the treated area to the sun within 48 hours after hair removal . Always pay attention to whether there is any bad skin reaction with sunlight exposure. Please use SPF30+ sunscreen for sun protection after 72 hours.
- Artificial skin tanning is not allowed during the hair removal operation.
- If you need to undergo plastic surgery after hair removal (only for the hair removal area), please consult your doctor.

Frequency for Reference

Time	Frequency	Fast Hair Removal Time	Fast Hair Removal Frequency
Week 1-4	Once / twice a week	Week 1-2	3 times a week
Week 5-8	Once a week	Week 3-4	Twice a week
After 8 weeks	Once half month / once a month	After 4 weeks	Once half month / once a month

Note: Adjust the use frequency according to the growth of the hair. For heavy hair, you can increase the level and frequency appropriately. The setting varies from person to person.

Cleaning and Maintenance

Please clean the product before and after each hair removal to avoid long-term accumulation of hair deposits and dust. If you do not clean normally, this product will not work as expected.

- Check this product, especially the light outlet port and body parts, before use. Make sure they are all clean and free from debris.
- Do not clean the product and its accessories in the sink or in the dishwasher.
- Do not use steel wool, abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, engine oil, acetone, etc.) to clean the product.

Note: If there are some spots you can no longer clean on the light outlet port, please stop using this product.

- Turn off the machine and allow it to cool before cleaning. If necessary, clean the surface of the product with a clean, soft cloth. Use the cloth to clean the light outlet port and ice compress panel.

Note: If the water cleaning effect is not good, you can use a few drops of high-concentration alcohol (not more than 96%) to clean the above parts.

Trouble Shooting

Problems	Possible Cause	Solutions
The device does not turn on	The adapter is not properly connected to the device or electricity supply.	Check the power connection
No flashing when in use	Temperature of main body is too high	After the device has cooled down, turn on it again
	The light outlet port is not completely fitted to the skin	
The display does not work		Contact consumer support

Problems	Possible Cause	Solutions
Feel painful when in use	Not shave the hair on the part to be treated	Before using this product, shave the hair on the part to be treated
	The part being treated is not suitable for using this product.	Do not use this product on the labia, vagina, anus, nipples and areola. Men should not use this product on the face or scrotum.
	Used a higher level	Choose to a lower level that you feel comfortable in use
	The light outlet port is broken	Contact consumer support
	Not use ice-cool mode	Press the ice-cool button to start the ice-cool mode
After the hair removal, the treated area turns red	It is normal that the treated skin may temporarily turns red	
The effect is not satisfactory	Used a lower level	Please choose a higher level next time
	Not flash on the area next to the previously treated area	Please flash in areas close to each other and there should be some overlap between flashing areas

Problems	Possible Cause	Solutions
The effect is not satisfactory	Insufficient frequency of using this device	To successfully remove hair, we recommend that you follow the reference frequency and level
	Your response to light exposure is slow	This product should be used continuously for at least 6 months as it can also inhibit hair regrowth during this period
The product becomes hot and humid during use	Just normal	No need to take any action
	The light outlet port is dirty or residual hair remained on the surface	Wipe the filter to remove the hair and reduce the light level

Error Code	Possible Cause
E1	Device overheating protection activated
E2	Head overheating protection activated
E3	The fan speed is too low or does not work
E4	The machine is broken
E5	The power adapter is broken

Warranty

Linglee warrants Hair Removal Machine U8 for 1 year from the date of purchase.

What's covered?

Defects due to faulty materials and/or workmanship of the product.

What's not covered?

1. Damage caused by the owner when attempting to fix or altering the product by himself or herself.
2. Damage caused by misuse, abuse, neglect, alterations or unauthorized repair.
3. Natural depreciation.

Warranty Period

Warranty becomes effective from the date of original purchase and will last for 1 year.

Who's covered?

This limited warranty is only effective for the original purchaser. It is non-transferable.

How to request warranty services?

To obtain the benefit of this warranty, please email linglee@mymemorism.com. Our customer service team will address your issue as soon as possible!

DISPOSAL OF USED PRODUCT BUILD-IN LITHIUM BATTERIES

Batteries may contain hazardous substances which could endanger the environment and human health if improperly disposed of. Batteries should be disposed of at an appropriate collection point for recycling. The recycling of materials will help to conserve natural resources. By ensuring the used batteries are disposed of correctly, you will help prevent potential negative consequences for the environment and human health. For more information about the recycling of used batteries, please contact your local municipality waste disposal service.



LINGLEE

www.mymemorism.com

linglee@mymemorism.com